

Low temperature plasma sterilizer for wide range of medical devices

The world's first direct injection pouch-type patented sterilization technology has been applied. SAL of 10⁻⁶ is validated to ensure sterilization stability of medical devices.

(SAL: Sterility Assurance Level, 99.9999%)



STERLINK Specification



Model	STERLINK FPS-15s Plus		STERLINK MINI	
Size (W x D x H)	433 x 614 x 437 mm		275 x 440 x 330 mm	
Chamber (W x D x H)	264 x 410 x 125 mm (14 L)		206 x 346 x 113 mm (7 L)	
Diagonal Length (Chamber)	470 mm		360 mm	
Vacuum Pump	Built-in type		Stand-alone type	
Weight	67kg		20kg (Pump module : 21kg)	
Cycle Time	Pouch Mode	STERPACK: 7 min	Pouch Mode	STERPACK: 7 min
	Pouch Plus Mode	STERPACK Plus: 14 min	Advanced Mode	STERLOAD MINI: 18 min
	Chamber Mode	STERLOAD: 36 min		

* Express mode is surface sterilization process, and other modes are penetration sterilization process.
 * All the sterilization process of STERLINK is double cycle.



sales@plasmapp.com
 www.plasmapp.com
 P-0001-2310-R0.9



STERLINK

7 Minutes Low Temperature Plasma Sterilizer



STERLINK

- Fast 7 min Sterilization
- Compact Design
- Easy Operation
- Real-time Monitoring



Smart Sterilization



Sterilant Cassettes, Consumables and Accessories

※ Sterilant cassettes



STERPACK
• Size: 135 x 280 mm
• 50 ea / Box



STERPACK Plus
• Size: 250 x 430 mm
• 30 ea / Box



STERLOAD
• 30 ea / Box



STERLOAD MINI
• 30 ea / Box

※ Consumables



Tyvek® 100
• Size: 100 x 400 mm
• 120 ea / Box



Tyvek® 200
• Size: 200 x 400 mm
• 90 ea / Box



Tyvek® 300
• Size: 300 x 400 mm
• 60 ea / Box



Sterilization Wrap
• Size: 60 x 60 cm
• 500 ea / Box

Accessories



STERSEAL™
• Size: 505 x 255 x 145 mm
• Weight: 12 kg
• Temperature: 135°C / 150°C



Cart
• Size: 483 x 660 x 603 mm
• Weight: 37 kg



Label Printer
• Size: 120 x 102 x 146 mm
• Weight: 0.5 kg



Label Sticker Roll
• Width: 52 x 110 mm
• Label counts: 49 ea

※ Accessories



CI tape
• Width: 20 mm
• Length: 50 m



CI strip
• Size: 18 x 105 mm
• 250 ea / Box



BI
• Time: 30 min
• 50 ea / Box

※ **Single Use Only**
※ Tyvek® is a Dupont™ registered trademark.

STERPACK & STERLOAD

Class IIb medical device certified by CE and TGA



Pouch Mode – 7 minutes fast sterilization with the world's first direct injection pouch and patented sterilization technology



Vacuum sealing pouches enable staff to have an immediate check on the sterile condition of the sterilized tools inside



Safe disposable sterilant cassettes

STERLINK FPS-15s Plus Specification

Model	Mode	Cycle time	Capacity	Sterilant
STERLINK FPS-15s Plus	Pouch Mode	STERPACK: 7 min	1 L	H ₂ O ₂ : 58%–59.5% (0.1 ml/cell)
	Pouch Plus Mode	STERPACK Plus: 14 min	4 L	H ₂ O ₂ : 58%–59.5% (0.3 ml/cell)
	Chamber Mode	STERLOAD: 36 min	14 L	H ₂ O ₂ : 58%–59.5% (0.9 ml/cell)

* STERPACK, STERPACK Plus and STERLOAD are certified sterilants by CE mark and TGA.

STERLINK™ MINI Specification

Model	Mode	Cycle time	Capacity	Sterilant
STERLINK™ MINI	Pouch Mode	STERPACK: 7 min	1 L	H ₂ O ₂ : 58%–59.5% (0.1 ml/cell)
	Express Mode	STERLOAD MINI: 12 min	7 L	H ₂ O ₂ : 58%–59.5% (0.7 ml/cell)
	Advanced Mode	STERLOAD MINI: 18 min		

* Express mode is surface sterilization process, and other modes are penetration sterilization process.
* All the sterilization processes of STERLINK are double cycles.

STERLINK verifies sterilization performance through the following lumen tests

Single-channel lumen claims for STERLINK

ø 0.7 x 500 mm Stainless steel | ø 2.0 x 1,500 mm Stainless steel | ø 1.0 x 2,000 mm PTFE

Compatible with most medical devices and materials that are heat and moisture intolerable



Simple to operate, install, and monitor (only requires electrical outlet)



Sterilization Performance

STERLINK verifies sterilization performance through the following lumen tests

Single-channel lumen claims for STERLINK


- ø 0.7 x 500 mm Stainless steel
- ø 2.0 x 1,500 mm Stainless steel
- ø 1.0 x 2,000 mm PTFE



Instrument Tracking System (ITS™)

Equipment tracking using barcode

- Real-time monitoring to check the STERLINK
- Provide the latest software updates remotely

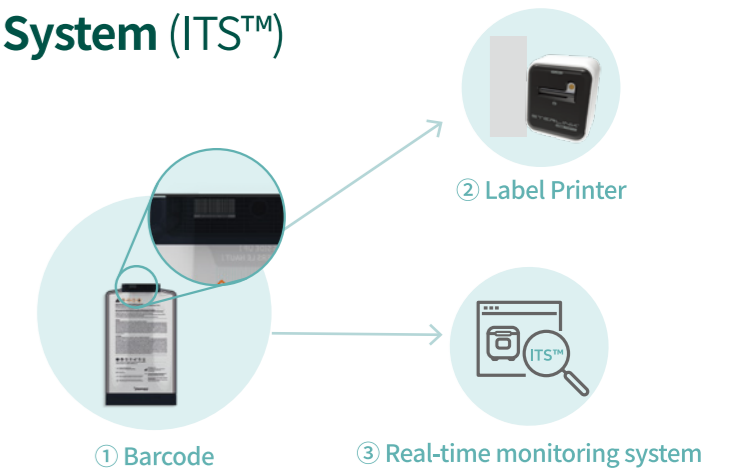

Comparison Chart

Sterilizer	Temperature	Cycle Time	Sterilant	Feature
Autoclave	Up to 134°C	60 min + 1 hour cooldown	Hot steam	• Cloth sterilization • Long cycle time • Risk of burn • Damage of medical instruments
E.O. Gas	Up to 60°C	Over 120 min + 12 hour ventilation	E.O. gas	• Highly dangerous toxic gas • Long cycle time
Plasma Sterilizer	Up to 60°C	70 min	H ₂ O ₂	• High cost • Large volume • Additional ventilation • Very low efficiency
STERLINK™	Up to 60°C	Pouch Mode : 7 min Pouch Plus Mode : 14 min Chamber Mode : 36 min	H ₂ O ₂ direct injection	• Fast sterilization cycle • Economic cost • Compact size • Eco-friendly • Ergonomic design • Easy maintenance • ITS™ system • Eco-friendly • Safeness of heat sensitive instrument